TAPIS DE MASSAGE EMS "SLIM"

FR

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FOOT MASSAGER

EN

Warning: Please read this instruction manual carefully before using the instrument and keep it for future use.

Warning:

The device works by the emission of electrical micro-pulses that are not necessarily perceptible from the start of using the product (depending on the degree of sensitivity of each).

It is advisable to choose intensity 1 from the start of the program, and, depending on how you feel after a few minutes, to increase the degree of intensity.

A. NOTES BEFORE USE

1. People meeting the following conditions should consult their doctor before using the product:
 - Patients extremely weakened by serious illness
 - For those with an integrated medical device, such as a pacemaker.
 - Suffering from heart disease, hypertension
 - Having recently undergone surgery.
 - For pregnant women or people who have recently been ill.
 - For people in need of medical treatment or who have had a previous physical abnormality being treated.
 - For people with injured feet or with an open wound.
2. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capacities, or by persons lacking in experience or knowledge, unless they are supervised or instructed on the use of this apparatus by a person responsible for their safety. Children should be supervised to ensure they do not play with this appliance.
3. Children should not use the product. They should be under adult supervision to ensure that they do not use the product as a toy.
4. Do not use on the face or on areas other than the feet.
5. Do not use while sleeping.
6. Do not allow this device to come into contact with jewelry or metal objects.

B. SAFETY INSTRUCTIONS


1. Only use the product as indicated in the user manual
2. Use the product only with the supplied accessories. Do not connect any other accessories.
3. Do not apply medical or cosmetic creams to the product or to your feet before use.
4. The recommended usage time is 15 minutes.
5. Make sure your skin is clear and clean before using the device.
6. Do not allow water to enter or splash the product.
7. Stop using immediately if you feel discomfort and nausea while using it.
8. Do not sit or stand on the device. Don't throw it away or drop it.
9. Contact of water with electricity can be dangerous. Do not leave the product in a humid environment (such as a bathroom or a swimming pool).
10. Do not use the product if the charging cable or the product itself is damaged.
11. In order to prevent possible electric shock, do not open the product.

12. If the device is damaged, do not attempt to repair it. This should be done by a specialist or technician. Otherwise, the guarantee will expire.
13. Do not use the product if the exterior of the product is damaged or if there is a mechanical or electronic failure.
14. Do not use the product in very hot conditions or near a heat source such as a heater.
15. Do not damage, bend or stretch the device which could cause possible electric shocks.
16. Be sure to use this device at a minimum distance of 1 meter from any device emitting microwaves, high frequency electromagnetism or short waves.

C. INTRODUCTION

Muscle electrical stimulation (EMS), also known as neuromuscular electrical stimulation, is induced by electrical impulses that produce multiple muscle contractions. The central nervous system contracts the muscles. EMS has proven to be an effective muscle building tool.

Interests of EMS: Muscle strengthening / Prevent muscle atrophy / Muscle relaxation / Relieve pain

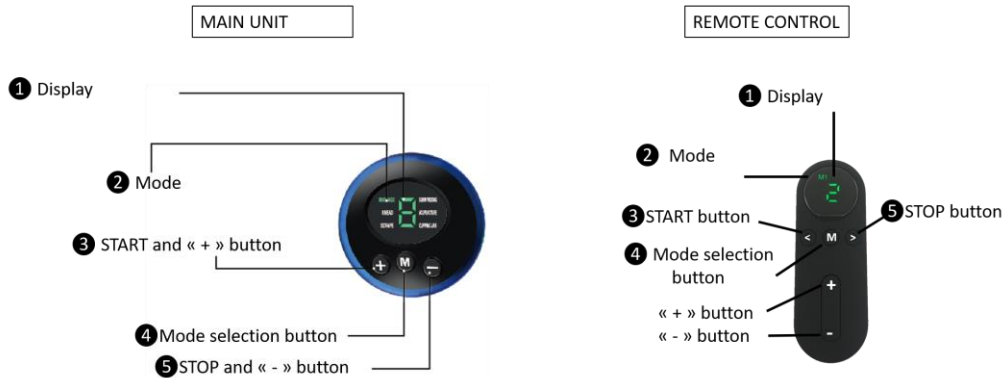D. HOW TO USE IT

1) Recharging the battery

The main unit and the remote control have a built-in battery. Before first use, charge them fully (charging time approximately 2 to 3 hours).

You can charge this device with a USB adapter (not included) or by connecting it to your computer's USB port (adapter specifications: 3V-1A)

The main unit indicator light turns red while the device is charging and turns off when the device is fully charged.

The remote will display "L" (low) while charging and then "F" (full) when charging is complete.

2) Start-up

This device has 6 massage modes and 9 intensity levels that can be adjusted via the main unit or the remote control.

- Clean your feet before using the appliance and dry them thoroughly and completely.
- Sit comfortably and place both feet on the mat surface.
- Start the device and select the appropriate mode and intensity setting via the main unit or the remote control. If you are using the main unit, after completing the settings you have 15 seconds to connect the main unit to the contact points of the foot mat.

If you are using the remote control, start the main unit, connect it to the contact points of the mat and then make settings with the remote control.

- Press the "ON" button **3** to switch on the device / remote control. The screen turns on.
- Use the "MODE" button **4** to select 1 of the 6 available modes among the following:


MAIN UNIT	Massage	Knead	Scrape	Elbow pressing	Acupuncture	Cupping jar
REMOTE CONTROL	1	2	3	4	5	6

- Then select the appropriate intensity by pressing the "+" button to increase the intensity and the "-" button to decrease it.

- The intensity of the stimulus is divided into 9 levels (0-9 lvl). When using for the first time, set the device to intensity 1 and then gradually increase according to your comfort.
- Enjoy a relaxing foot massage.

This appliance is automatically set to a 15 minute timer, so it will automatically turn off when the time has elapsed.

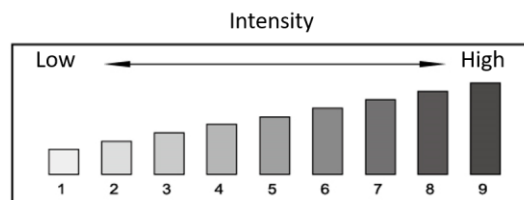E. CARE AND MAINTENANCE

Clean the mat surface after each use.

Use a damp cloth or cloth to wipe the surface of the carpet.

Do not rinse or immerse the device in water.

Do not use abrasive cleaners or solvents.

If you will not be using this device for an extended period of time, store it safely in a clean, dry place.

Do not disassemble or destroy this device.

F. FEATURES

- Materials: PU

- Power supply: 240 mah lithium battery and 150 mah rechargeable battery for the remote control included (1 USB charging cable included). DC 3V - Frequencies 1 ~ 100Hz

- Session: 15 minutes

- Autonomy: 2 h (about 8 cycles of 15 min)

- Charging time: 2-3h

The European directive 2012/19/EC about electronic and electric waste, requires that you can't throw away defective domestic appliance with common waste. Used device has to be collect separately to optimize the recuperation rate and the recycling of the materials to reduce the impact on the health and the environment.

Precision and precautions about the battery use

- *Throwing away a used battery in the environment or with household waste pollutes and avoids the recuperation of reusable materials.*
- *It is important to limit the battery consumption and to respect the following rules:*
 - *Favor alkaline battery (last longer than saline battery) and when it's possible use rechargeable battery.*
 - *Leave used batteries and accumulators in adapted containers that you can find in some retailers. In that way, metals don't pollute the environment and the health.*
- *Battery has to be inserted by respecting the polarity show on the apparel and the battery. An incorrect positioning can damage the apparel, cause leaking and even provoke fire and the explosion of the battery.*
- *To assure a good functioning, battery has to be in a good condition. If there is a problem, put new battery.*
- *Never try to recharge battery not made for this. It could leak, warm up, cause fire or explode.*
- *Replace all the battery at the same time. Never mix up alkaline battery with saline battery.*
- *Used battery has to be removed.*
- *Also, remove the battery if you're not using the apparel for a long time.*
- *Never try to short-circuit the terminal of a battery.*
- *Never throw away the battery in fire because they can explode.*
- *Take out the accumulator from the apparel before reload it.*
- *We recommend that an adult watch if a child change the battery to respect those rules.*
- *If you swallow a battery, go immediately to your doctor or the hospital.*
- *For the best ongoing performance, we recommend using only fresh, high-quality, alkaline batteries.*

Importé par / Importiert von / Geïmporteerd door / Imported by :
CAPTELEC, 59170 Croix - FRANCE

Carton et emballage papier à
 trier / Kartons und
 papierverpackung recyceln /
 Verpakkingsmateriaal in karton
 en papier te sorteren / Recycle
 carton and paper packaging